



## PCP01, PCP02 & PCP03 Vest and Leotard Order Form

All fields are required in order to process your order

### Order Details

Date: \_\_\_\_\_ Date Required: \_\_\_\_\_  
Order No.: \_\_\_\_\_  
Contact Name: \_\_\_\_\_  
Contact Phone No.: \_\_\_\_\_  
Email: \_\_\_\_\_  
Hospital/Clinic: \_\_\_\_\_  
Delivery Address: \_\_\_\_\_  
\_\_\_\_\_  
\_\_\_\_\_  
Post Code: \_\_\_\_\_

### Patient Details

Patient Reference No.: \_\_\_\_\_  
Is this their first SDO®?  No  Yes  
First Name: \_\_\_\_\_  
Surname: \_\_\_\_\_  
Date of Birth: \_\_\_\_\_  
Diagnosis: \_\_\_\_\_  
\_\_\_\_\_  
Contact me about my order via:  Email  Phone  
Attached with order:  Images  Consent form\*  
\*(Signed consent form must be supplied when providing patient images)

Clinician Number: \_\_\_\_\_

### Specifications

Garment:  Vest  Leotard  
Sleeves:  No sleeves  Short sleeves  Long sleeves  
Crotch: (Leotard only)  Poppers  Velcro  Standard Fit  Higher (see drawing on page 5)  
Gastrostomy site:  No  Yes (Please mark position at first fitting only, if a gastro opening is required, the garment will be supplied without front zip)  
Silicone Edging:  No  Yes  
Velcro tabs:  No  Yes (To attach to leggings)  
Soft Binding:  Yes (Soft binding will be applied to neck and crotch only, in complementing colour or theme if selected)

### Zip Position (please specify)

Front:  No zip  Closed  Open  Upside down  Reflux / Abdominal Zip  
Back:  No zip  Closed  Open  Upside down  
Sleeves:  No zip  Full length  Back  Mid forearm

### Themes

1 Spaceman  2 Astronaut  3 Footballer  4 Striker  5 Goalie  6 Noir  
 7 Silver Aztec  8 Pink Aztec  9 Leopard Flower  10 Lavender  11 Daisy  12 Ultraviolet  
 13 Mermaid  14 Starlight  15 Marshmallow  16 Tribal  17 Blue Frenchie  18 Fawn Frenchie  
 19 Cadet  20 Lance Corporal  21 Pilot Officer  22 Admiral  23 Major  24 General

Please select one of the themes above, your garment will be manufactured according to the theme selected. Alternatively, select a plain coloured garment below. Without either a theme or plain garment selection, we cannot proceed with your order.



### Plain Coloured Garments (Base fabric, panels, binding and thread are all the same colour)

White  Beige  Pink  Black  Navy Blue/Black\*  
\*Navy base fabric, black panels, binding & thread

When completed, please click: [customerservice@jobskin.co.uk](mailto:customerservice@jobskin.co.uk) to email your electronic order form

Please download your electronic forms directly from our website - [www.jobskin.co.uk/file-download](http://www.jobskin.co.uk/file-download)

Order No.: \_\_\_\_\_ Patient Reference No.: \_\_\_\_\_

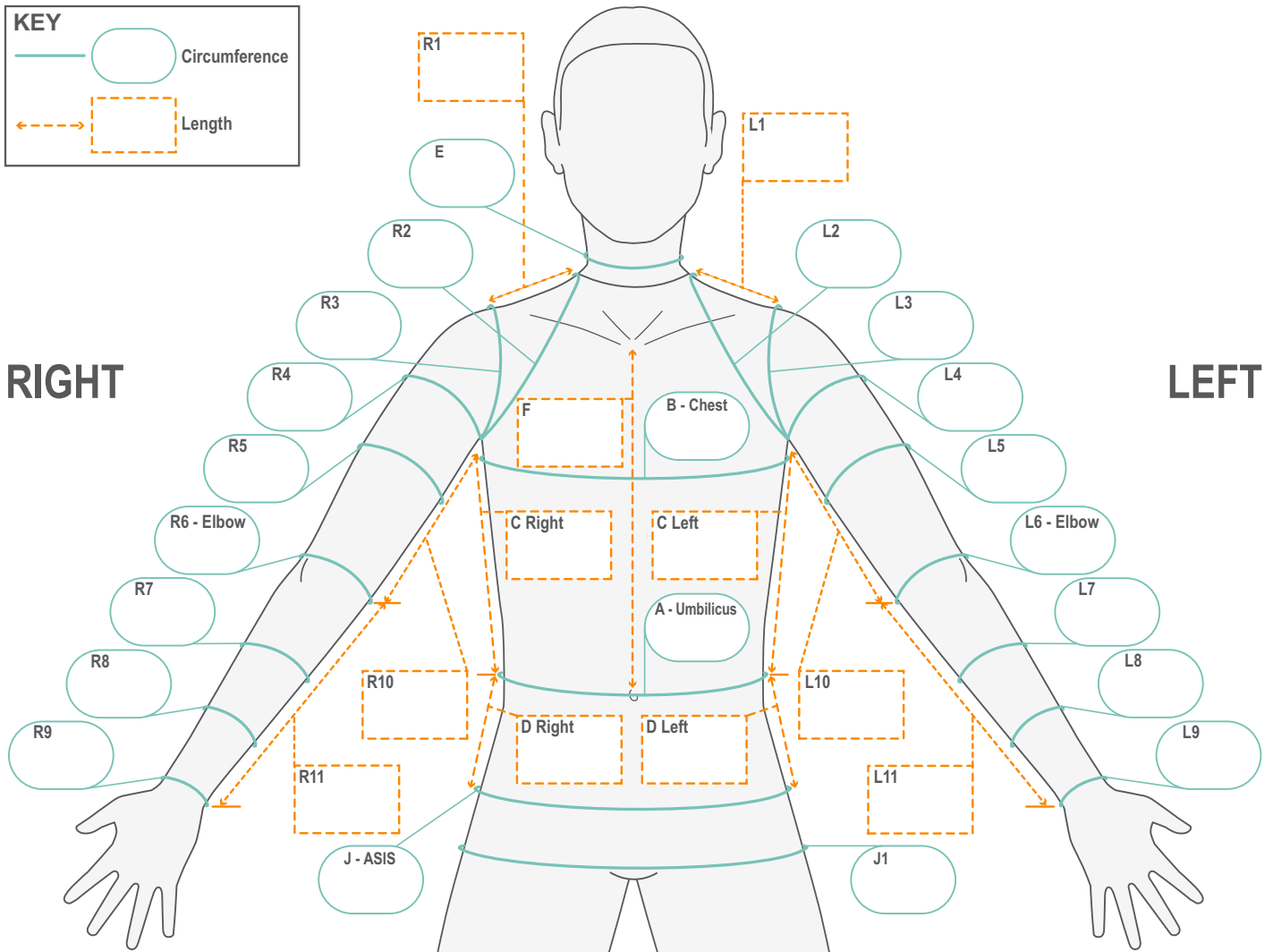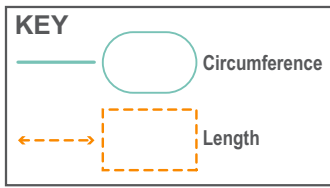

**Important:** Please make sure when measuring for an SDO® garment that the patient is measured lying down, supine and in the best position of alignment.

### Circumference Measurements (C)



			Left (cm)	Right (cm)
A	C	Level with umbilicus, 10mm less than 'at rest' circumference		
B	C	Level with chest in line with the axillary crease		
E	C	Base of the neck		
J	C	Level with Anterior Superior Iliac Spine (ASIS) or end of garment		
J1	C	Level with widest part of buttocks		
2	C	Shoulder joint from base of the neck (at point where neck meets shoulder below ear lobe), under the axilla and back to the same point. (Scapula set against the chest wall and the arm by the side of the body). One should be able to slide a finger under the tape measure at the front of the axilla		
3	C	Shoulder joint acromion process under the axilla and back to the same point. (Scapula set against the chest wall and the arm by the side of the body). One should be able to slide a finger under the tape measure at the front of the axilla		
4	C	Upper arm level with axilla. Arm should be by the side		
5	C	Mid upper arm level with muscle bulk of biceps, 5mm less than 'at rest' circumference		
6	C	Elbow joint with arm extended or end of sleeve		
7	C	Forearm level where muscle bulk is greatest, 5mm less than 'at rest' circumference		
8	C	Forearm at musculotendinous junction (approximately 1/3 of forearm up from the wrist)		
9	C	Wrist level with wrist crease or end of sleeves		

Order No.: \_\_\_\_\_ Patient Reference No.: \_\_\_\_\_

### Length Measurements (L)



Left (cm)

Right (cm)

C	L	Point level with axillary crease measure down side to a point level with umbilicus		
D	L	Waist to hips or end of garment		
F	L	Below sternal notch to umbilicus on the anterior surface of the chest with the ribs held in best position of alignment. This measurement dictates the scoop of the neck at the front		
10	L	Axilla to cubital crease on the medial border or to end of short sleeve		
11	L	Point level with the cubital crease to wrist crease or end of long sleeve along ulnar border of forearm		

### Additional Measurements Lengths (L)



Left (cm)

Right (cm)







1	L	Acromion to base of neck or width required (only for Vest or Leotard with no sleeves)		
---	---	---	--	--

### Reinforcement Panels (Please specify)

<p><b>Standard Panels</b></p> <p>To assist hip external rotation and back extension</p> <p><b>SP</b></p> <p>Std panels included with Suits. Type NO if no panels are required _____</p> <p>Posterior</p>	<p><b>Lumbar Panels</b></p> <p>To assist postural control</p> <p><input type="checkbox"/> <b>LP</b></p> <p>Extend to _____ cms above umbilicus</p> <p>Anterior    Posterior</p>	<p><b>Lateral Trunk Panels</b></p> <p>To correct left lateral flexion</p> <p><input type="checkbox"/> <b>TPL</b></p> <p>Anterior    Posterior</p>	<p><b>Lateral Trunk Panels</b></p> <p>To correct right lateral flexion</p> <p><input type="checkbox"/> <b>TPR</b></p> <p>Anterior    Posterior</p>	<p><b>Chest Panels</b></p> <p>To assist protraction of the scapula</p> <p><input type="checkbox"/> <b>CP1</b></p> <p>Anterior</p>
<p><b>Cross Panels</b></p> <p>To assist retraction of the scapulae and back extension</p> <p><input type="checkbox"/> <b>CP2</b></p> <p>Posterior</p>	<p><b>Abdominal Panels</b></p> <p>To resist lordosis</p> <p><input type="checkbox"/> <b>AP</b></p> <p>Anterior</p>	<p><b>Back Panels</b></p> <p>To resist kyphosis</p> <p><input type="checkbox"/> <b>BP</b></p> <p>Posterior</p>	<p><b>Oblique Panels</b></p> <p>To assist activation of oblique muscles</p> <p><input type="checkbox"/> <b>OP</b></p> <p>Anterior    Posterior</p>	

Order No.: \_\_\_\_\_ Patient Reference No.: \_\_\_\_\_

## Scoliosis Reinforcement Panels (Please specify)
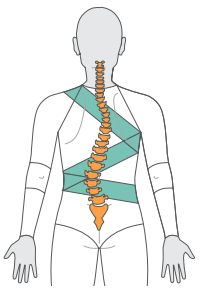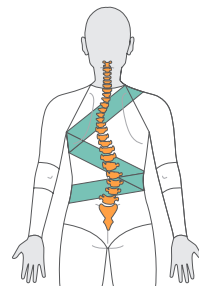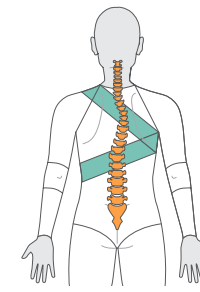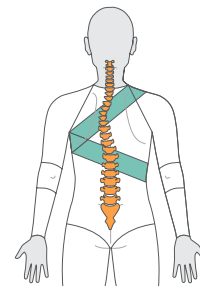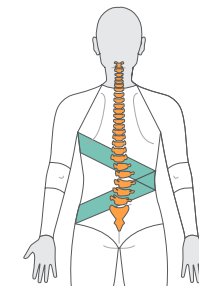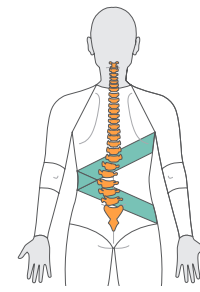
<p><b>S-Shape Scoliosis Panels</b></p> <p>To assist correction of right thoracic and left lumbar curves</p> <p><input type="checkbox"/> SPL</p>  <p>Posterior</p>	<p><b>S-Shape Scoliosis Panels</b></p> <p>To assist correction of left thoracic and right lumbar curves</p> <p><input type="checkbox"/> SPR</p>  <p>Posterior</p>	<p><b>C-Shape Scoliosis Panels</b></p> <p>To assist correction of right thoracic curve</p> <p><input type="checkbox"/> CPL</p>  <p>Posterior</p>	<p><b>C-Shape Scoliosis Panels</b></p> <p>To assist correction of left thoracic curve</p> <p><input type="checkbox"/> CPR</p>  <p>Posterior</p>	<p><b>C-Shape Scoliosis Panels</b></p> <p>To assist correction of right lumbar curve</p> <p><input type="checkbox"/> Low CPL</p>  <p>Posterior</p>	<p><b>C-Shape Scoliosis Panels</b></p> <p>To assist correction of left lumbar curve</p> <p><input type="checkbox"/> Low CPR</p>  <p>Posterior</p>
---	--	---	---	---	--


### Optional Scoliosis Measurements Left (cm) Right (cm)

L	Drop from axilla down to centre point of first panel		
L	Waist up to centre point of second panel		
L	Drop from hip to centre point of last panel		

### Optional Scoliosis Measurements (cm)

L	Width of panel *	
---	------------------	--

\* NOTE: Paediatrics - width of panel must be between 4-10cm  
Adults - width of panel must be between 6-12cm

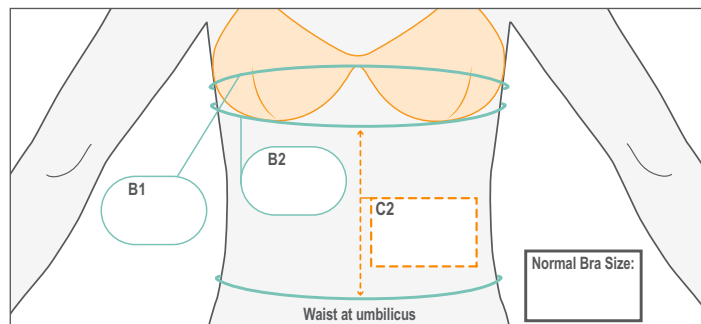


### Measurements for Bra Vest (cm)

B1	C	Overbust circumference	
B2	C	Underbust circumference	

C2	L	Length from underbust to waist	
----	---	--------------------------------	--

Normal bra size		
-----------------	--	--

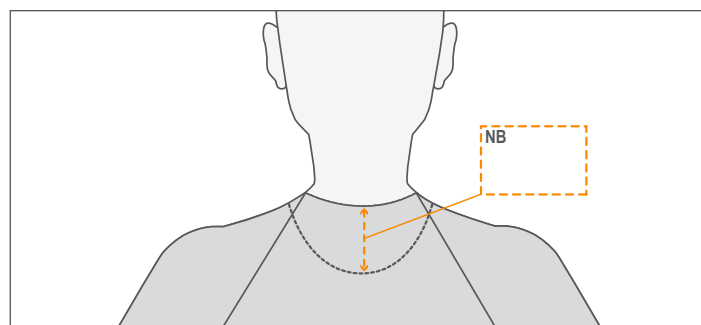


### Back Scoop Neckline Options (cm)

NB	L	Neckline scoop length at back	
----	---	-------------------------------	--

For neckline scoop at front, see F measurement

### Style Options

V Neck       Round Neck

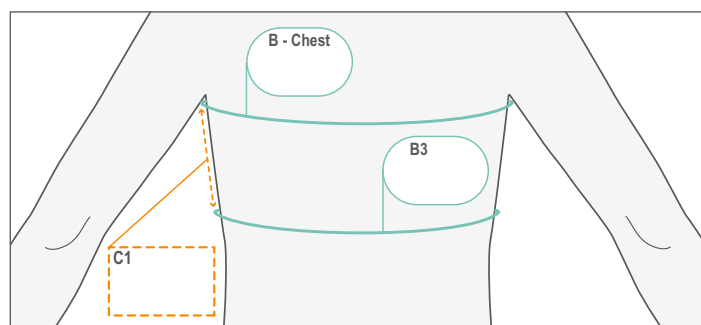


### Measurements for Rib Flare (cm)

B	C	Circumference level with chest in line with the axillary crease	
B3	C	Rib flare circumference at widest part	


C1	L	Length from axilla to rib flare at widest part	
----	---	--	--

Order No.: \_\_\_\_\_ Patient Reference No.: \_\_\_\_\_

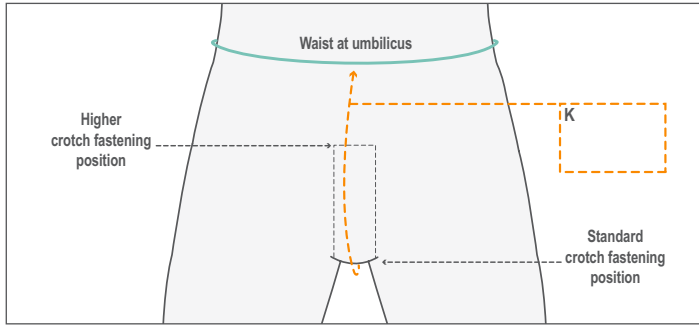

### Measurements for Leotard (cm)

K	L	Waist to waist through legs	
---	---	-----------------------------	--

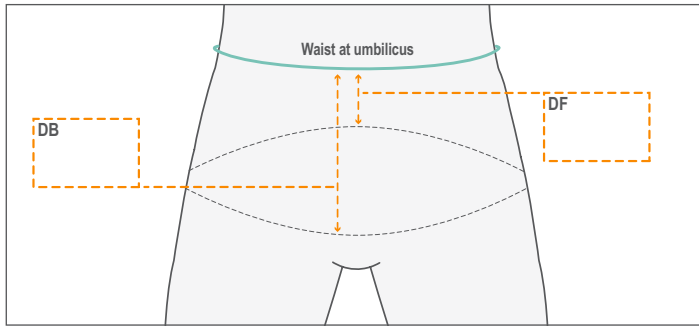

### Measurements for Back & Front Scoop (cm)

DF	L	Waist down to end of garment at the front	
DB	L	Waist down to end of garment at the back	